

Vet-PDA Occluder™

Implantable device for closing a persistent arterial duct



Vet-PDA Occluder™

The Vet-PDA Occluder™ is a permanently implantable device designed to safely and effectively close a persistent arterial duct. Manufactured with high-quality biocompatible materials, the Vet-PDA Occluder™ is specifically designed for the anatomy of small animals such as dogs and cats weighing less than 3 kg, ensuring proper function and minimizing potential complications.

Benefits:



- **Efficiency:** The device ensures a secure and permanent closure of the persistent arterial duct, improving blood circulation and the animal's health.

- **Safety:** Made with biocompatible materials like Nitinol, the device minimizes the risks of rejection or postoperative complications.



- **Veterinary Adaptation:** Enhanced coil flexibility, same 4Fr introducer catheter with an optimized length of 55cm, and a curved tip for better alignment during implantation maneuvers.

- **Ease of Implantation:** Its design allows for simple and quick placement, reducing surgical time and associated risks through a non-invasive approach.




- **Faster Recovery:** Due to its easy placement and minimal invasion, the device enables quicker and less painful recovery for the pet, reducing hospitalization time and potential postoperative complications.

Technical Features:

- Nitinol spiral disc for arterial ducts up to 4 mm.
- Pre-assembled system.
- Distal reinforcement to prevent recoil.
- Release system with radiopaque markers.
- Flexible and adaptive design for precise and safe coil release.
- Repositionable and retrievable before release.
- Immediate closure.
- Single-use.
- Suitable for all types of PDA.

Vet-PDA Occluder™

Data Sheet:

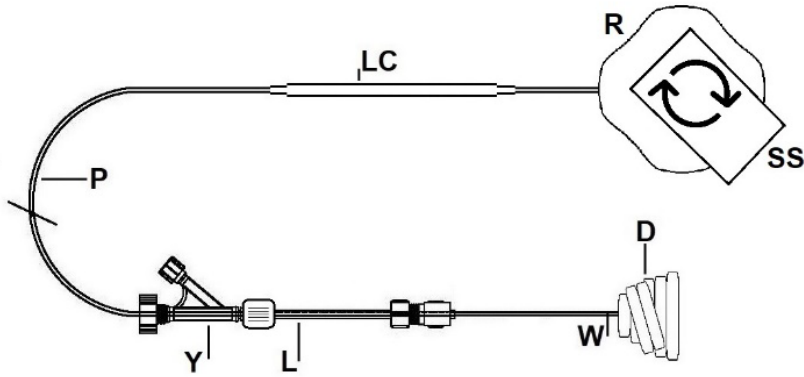

Tradename	Vet-PDA Occluder™
(D) Device	Nitinol SE 508
(W) Guide Wire	Nitinol SE 508
(L) Loader	Polycarbonate, HDPE, PTFE
(Y) Y Connector	Transparent Polycarbonate
(P) Pusher	Stainless Steel 1.4301
(LC) Loading Catheter	HDPE w/ 20% BaSO4
(R) Release Handle	ABS
(SS) Security Seal	Self-adhesive Paper
Angiographic Catheter	Pebax
Sterilization	Ethylene Oxide (EO)
Packaging	Individual Sterile Packaging
Conservation	Normal Temperature and Humidity
Manufacturer	Evomed S.L.U. Ctra. Torrejón-Ajalvir Km. 5,2 28864 Ajalvir (Madrid) Spain

Vet-PDA Occluder™

References:

Reference	Distal x Prox. Coil Diameter	Type	Coil Length (mm)	Implantation Catheter	Catéter Length
VETPDA000404	4 x 4 mm	FLEX	3.5	4 Fr	55 cm
VETPDA000504	5 x 4 mm	FLEX	3.5	4 Fr	55 cm
VETPDA000605	6 x 5 mm	FLEX	3.5	4 Fr	55 cm
VETPDA000706	7 x 6 mm	FLEX	4.5	4 Fr	55 cm

The Vet-PDA Occluder™ is the preferred choice of veterinarians worldwide due to its proven effectiveness and exceptional results in treating persistent arterial ducts in dogs and cats. With cutting-edge technology and high-quality materials, the Vet-PDA Occluder™ ensures the best solution for securely and effectively closing the anomalous duct, providing pets with a healthier and more active life.

Contact:

Distributor

Evomed S.L.U.
Carretera Torrejón-Ajalvir Km. 5,2
28864 Ajalvir, Madrid
Spain
Tel: +34 918 874 220
www.evomedvc.com
info@evomedvc.com

Manufacturer

Evomed S.L.U.
Carretera Torrejón-Ajalvir Km. 5,2
28864 Ajalvir, Madrid
Spain
Tel: +34 918 874 220
www.evomedvc.com
info@evomedvc.com